



**RELIABLE – QUIET – EASY TO USE**

# Quattro

## Bedpan Disposal System

---

*With proven reliability, easy to clean surfaces, class-leading quietness and efficiency, the **Quattro** by Haigh has been designed with Nurses, Estate Managers and Site Operators front of mind*

---



## Reliability

We absolutely understand that your bedpan machine must work without fail tens, if not hundreds of times per day in a busy sluice-room.

**Quattro** is tested on pulp disposable items from all of the many manufacturers and types of product, from urine bottles to large wash-bowls and sample pots

Haigh test all machines rigorously as they leave our production line. Many features have been expertly built in specifically to reduce the risk of machine failure or blockage of your drainage system:

- ✓ Only particles < 5mm x 8mm can exit hopper
- ✓ Automatic and mechanical operation
- ✓ Ease of access on all key components

Foreign objects (gloves etc.) are held back in the **Quattro**, preventing blockages in your site drainage.

**FAST, QUIET, CLEAN, EASY**

**SUPERFINE MACERATION**

**PROTECT YOUR DRAINS**

**QUALITY ASSURED ON ALL TYPES OF PULP**


## Compliance

Your sluice room may not be the most glamorous part of your site but it need not be unpleasant to work in and needs to enable compliant with relevant regulations.

**Quattro** improves your site environment:

- ✓ Easy to clean
- ✓ Embedded antimicrobial silver ions
- ✓ Efficient use of Energy and Water for ISO14001
- ✓ Hopper seal contains odour and aerosols
- ✓ End of cycle sanitization
- ✓ Anti-legionella system
- ✓ WRAS approved

Cleaner surfaces combined with the single use pulp concept demonstrably leads to lower risk of cross infection.

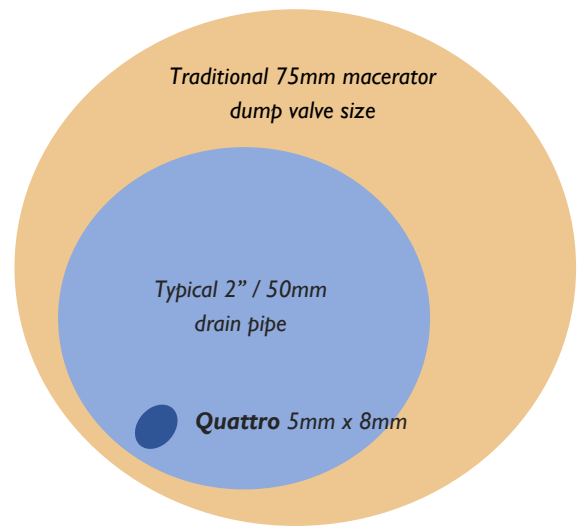

## Nice to Use

**Quattro** is designed to make the life of the user as straightforward and pleasant as possible. As the 5<sup>th</sup> generation of the Haigh Sluicemaster range **Quattro** does this by giving you:

- ✓ The lowest loading height
- ✓ The widest opening
- ✓ High visibility foot pedal opening
- ✓ Automatic starting
- ✓ Fast cycle time for the 'morning rush'
- ✓ The quietest macerator as tested by Quietmark

As a result, not only does Quattro aim make the life of the carer better but particularly with its class-leading low sound levels and short run time, the disposal of patient waste is kept discrete with minimal disruption, even during the night.




## **NATIONWIDE SUPPORT**

*Ask us about:*

24hr Response | Maintenance Agreements  
Genuine Spare Parts | Vernacare Maintenance  
Annual Contracts | Extended Warranties  
Installations | Asset surveys | Commissioning  
Decommissioning | Replacement Schedules

**01989 763131    [info@haigh.co.uk](mailto:info@haigh.co.uk)**

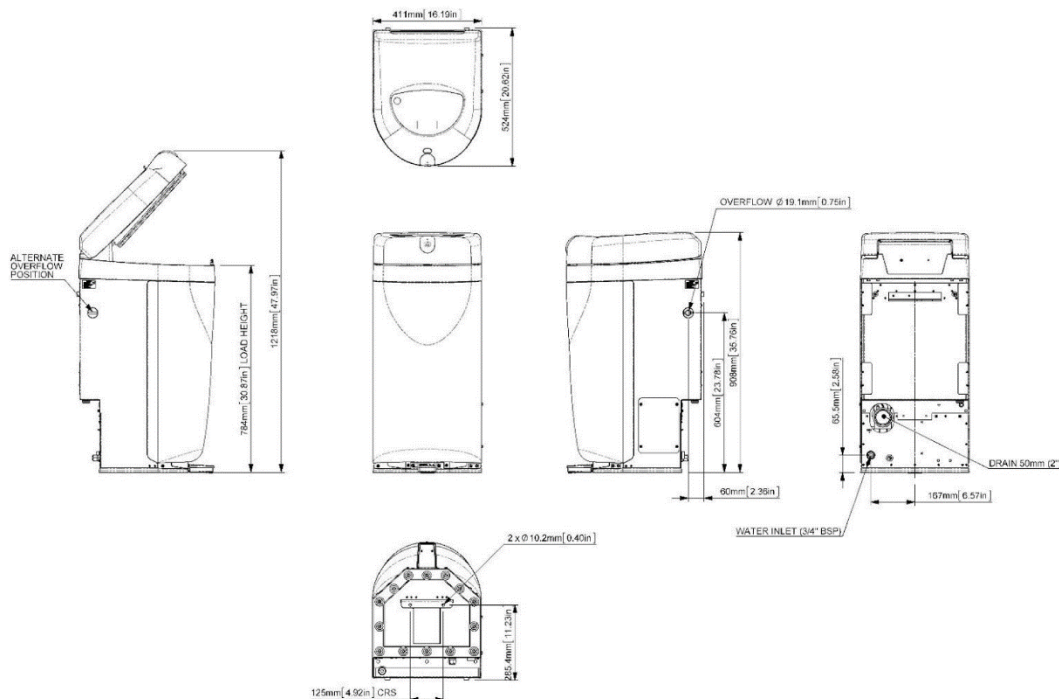

# FAST – QUIET – EFFICIENT

## Technical Performance

Cycle Characteristics: 90-180 seconds | 52-56 dBA | 54-59 dB max | Quietmark approved | Touch Free Open

Water: 15 - 22 litres per cycle used | 5.5 litres per minute cold water supply required

Electricity: < 0.01 kWh per cycle | 230V / 1Ph / 50Hz or 400V / 3Ph / 50Hz supply required

Inlet, Outlet & Cistern: 3/4" BSP Inlet | 2" BSP multi-fit P-trap Outlet | 11 litre cistern: EN 13077

Metrics: Width: 411 mm | Loading Height: 784mm | Depth: 524mm | Dry Weight: 58.5 kg

Haigh Engineering Ltd,  
Alton Road,  
Ross-on-Wye,  
Herefordshire, HR9 5NG,  
United Kingdom

Web: [www.haigh.co.uk](http://www.haigh.co.uk)

Tel: 01989 763131

Email: [info@haigh.co.uk](mailto:info@haigh.co.uk)

